

Changing Places Product Schedule

Changing Places Overview

CHS Healthcare is a leader in the health care industry, providing quality and innovative patient transfer solutions and hygiene systems to the Australian market. At CHS Healthcare we go above and beyond to ensure that the best outcomes are achieved in a simple, concise and cost effective way.

We understand the market and have produced this schedule to help designers, planners and the construction industry navigate the complex Changing Places design and construct process. CHS Healthcare is proud to have been involved in the Changing Places initiative since its introduction to Australia in 2012. We are actively involved on a national level, supplying and installing equipment in every state of Australia. Our expertise in this field is valued by clients and Changing Places alike, resulting in CHS Healthcare being crowned 'Industry Champions'.

We look forward to helping you out on your Changing Places project in anyway possible. Even if we don't have access to what you need we have key relationships with numerous associated suppliers, architects, access consultants and other valuable contacts to make your job simpler.

Kind Regards,

Andre Jensen | Project Consultant
P: 1300 789 420 | M: 0448 741 367
E: andre@chshealthcare.com.au | W: www.chshealthcare.com.au

Signature: _____




Changing Places Hoist Ceiling Clearances


A. SUPA RAIL PARALLEL RAIL AND SUPA RAIL TRAVERSE RAIL - XY TROLLEY MANUAL STANDARD



B. SUPA RAIL PARALLEL RAIL AND SUPA RAIL TRAVERSE RAIL - XY TROLLEY ABOVE MOUNT

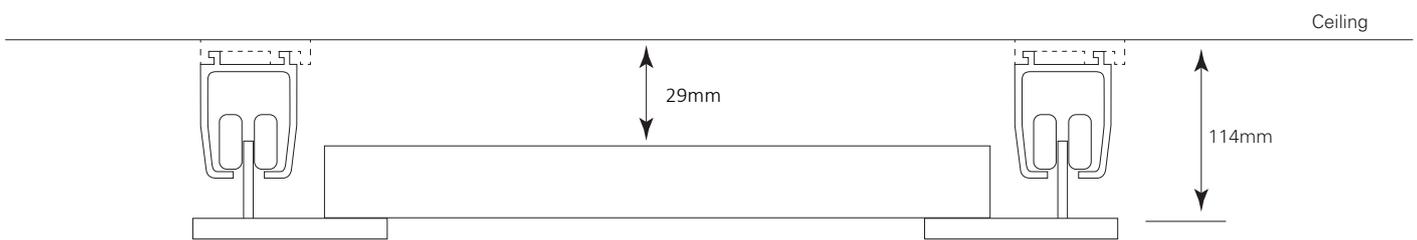

C. FINELINE (CONCEALED) PARALLEL RAIL AND SUPA RAIL TRAVERSE RAIL - XY TROLLEY MANUAL STANDARD

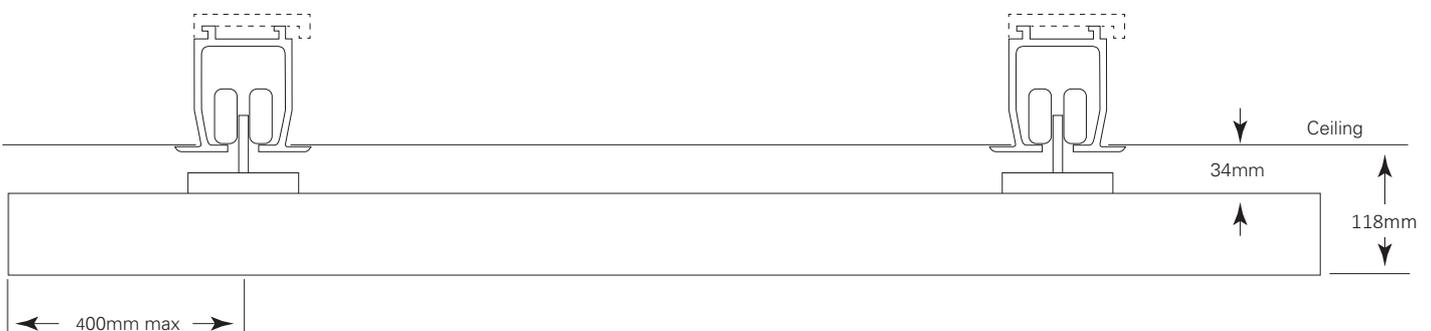

Need more information?


web

email

call

SUPARAIL

Rail Systems

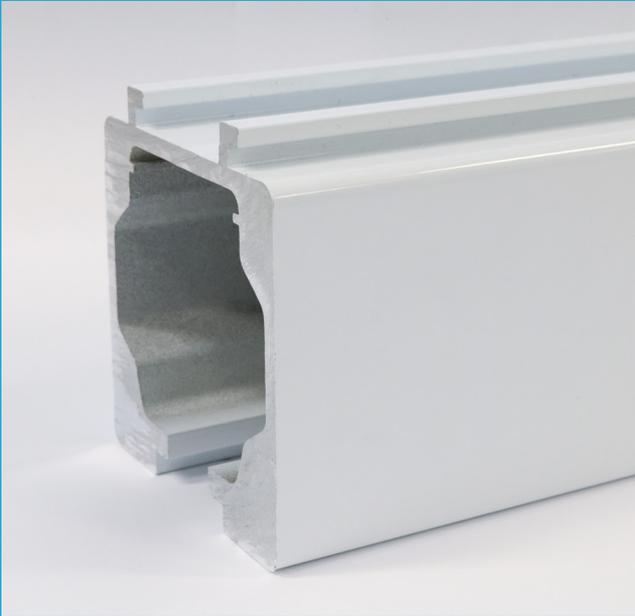

The thicker profile of the SUPARAIL allows for extended spans between supports making it suitable for XY traverse rails and wall mount situations. Reduced depth provides a more aesthetically pleasing finish without loss of lift height.

Technical Specifications


- Height 84mm
- Width 76mm

Product Code:

3.0M - CHSSR03
4.0M - CHSSR04
6.0M - CHSSR06

Features

- Suitable for extended spans
- Suitable for XY traverse & wall mount situations
- Reduced depth allows for polished finish

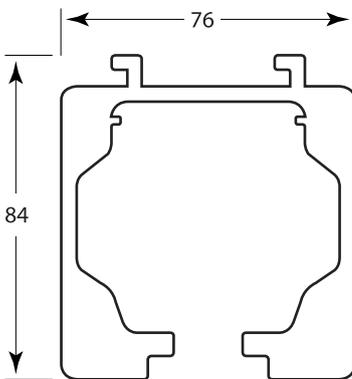

Need more information?

web

email

call

CHS Fineline Rail

Rail Systems



CHS Fineline rails are designed to be built into closed ceilings, giving a beautiful integral appearance. This concealed lightweight rail is suitable for systems up to 450kg SWL and suits CHS's full range of fixed and portable hoist units. Neat moulded access panels at rail end provide for hoist servicing needs

Technical Specifications

- Height 76mm
- Width 58mm
- Track Opening 18mm
- Base Width 81mm

Product Code:

2.0M	3.0M	4.0M	6.0M	7.2M
LFL02	LFL03	LFL04	LFL06	LFL072

Features

- Suitable for systems up to 450kg SWL
- Suits all CHS Fixed & Portable Hoist Range
- Concealed appearance

Need more information?

 web email call

440F Fixed Hoist



200kg SWL Fixed Hoist

The CHS Stellar Series ceiling hoist offers a cost effective, flexible solution. The hoist has the new Kinetic Charge System, and automatically recharges whenever the hoist is lowered. The 440F is equipped with a quick and safe emergency lowering system. SWL of 200Kg. On-rail service capability. Powered traverse is also available.

Technical Specifications

- 24V Motor
- Emergency Stop & Emergency Lowering
- 2 x 12V Batteries
- In-built service counter & diagnostics
- Lift Speed 42mm/sec

Features

- Hybrid automatic recharging system allows hoist units to automatically recharge, giving more battery lift than any other ceiling hoists
- Heavy duty capacity system with single lifting motor, allows for lifting, transferring & positioning of virtually every patient
- Manual or power traversing models available
- Strong Build Quality
- Quick change battery pack whilst on the rail
Constant charge options available for Changing Places applications

Need more information?

 web email call

Waverley Glen C450

204kg SWL Fixed Hoist



The Waverley Glen C450 Hoist delivers innovation, affordability and quality. Ceiling hoists dramatically reduce injuries resulting from the handling of patients and residents in both institutional and home care settings. The C450 has been designed with a full set of standard features and can be safely operated by a single caregiver. Built for quiet operation, the C Series provides a smooth & effortless transfer. The C Series will lift up to a maximum of 205kg (450lbs), allowing greatest flexibility in the type of patients that can be transferred. For additional flexibility, powered traversing options can be included.

Technical Specifications

- Lift Motor 24V DC
- Traverse Motor 24V DC
- Charger Input 240V AC, 1.0 Amps
- Charger Output 24V DC, 1.5 Amps
- Battery 24V (2 x 12V DC), 5.0AH
- Hoist Case Flame Retardant ABS
- Hand Control Pneumatic 100% waterproof
- Yoke Stainless Steel
- Lifting Height up to 213cm
- Unit Weight Under 10kgs
- Max Load 204kg
- Lift Speed 46mm/sec

Product Code: 323100 / 323149

Features

- Smallest Hoist of its class which allows for 180mm of additional lifting height
- Quick lifting speed
- Emergency pull cord.
- Modular design for simplified maintenance
- Quick charge batteries for reduced service time
- Audible low battery indicator & visual battery/charge level display
- Digital readout
- Auto shut of when not in use to save battery life
- 2.2 strap vertical movement allows you to easily lift from the floor
- Pneumatic Hand Controls – 100% waterproof
- Manual option – for emergency raising or lowering
- Available with Return to charge option.

Need more information?

 web email call

Nivano CTE Changing Table

Change Tables




Our easy-to-operate multi-purpose height-adjustable shower and changing tables come in standard lengths from 1200 to 1900mm, with a weight capacity of 200kg. Non-standard sizes are also available.

The Nivano changing table has a travel range of 700mm and can easily lift a client up to 200kg from near floor level to a convenient height for washing and changing. It folds away very neatly when not in use, protruding just 220mm from the wall and features battery back-up for safe lowering in the event of a power failure.

Technical Specifications


- Standard Lengths from 1200 – 1900mm
- Weight Capacity 200kg
- Travel Range 700mm
- Folded Width 220mm

Product Code: CTE

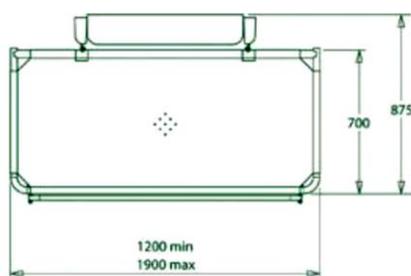

Features

- Folds neatly away
- Battery back-up for safe lowering
- Lifts from near floor level
- Optional fold down guard and guard bumper
- Optional removable mattress
- Optional padded top
- Optional anti-splash skirt
- Optional tilting backrest

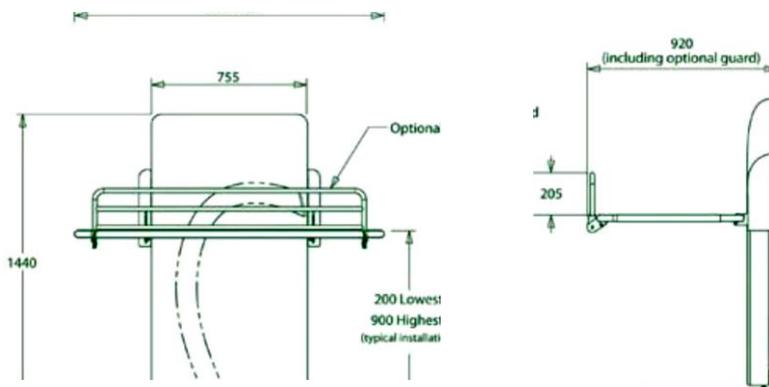

Need more information?


web

email

call

Nivano CTX Changing Table

Change Tables

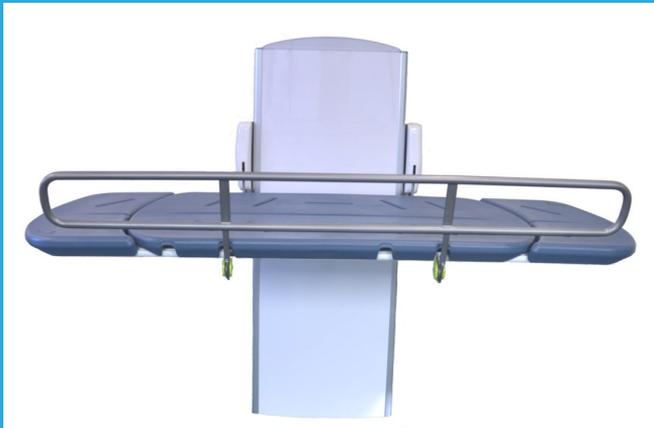

Our easy-to-operate multi-purpose height-adjustable shower and changing tables come in two standard lengths of 1500 and 1800mm, with a weight capacity of 200kg. The Nivano CTX changing table has a travel range of 700mm and can easily lift a client up to 200kg from near floor level to a convenient height for washing and changing. It folds away very neatly when not in use, protruding just 310mm from the wall and features battery back-up for safe lowering in the event of a power failure.

Technical Specifications

- Standard Lengths 1500 & 1800mm
- Weight Capacity 200kg
- Travel Range 700mm
- Folded Width 310mm

Product Code: CTX

Features

- Folds neatly away
- Battery back-up for safe lowering
- Lifts from near floor level
- Optional fold down guard and guard bumper
- Integrated water collection tray
- Moulded Polyurethane padded top
- Integrated tilting backrest/headrest

Need more information?

 web email call

CHS Change Table Cover Dispenser



CHS adult change table cover dispenser is designed to provide a clean and sanitary way to protect the surface of a change table. Meeting all AUS standards and the Changing Places kit guidelines this is perfect any application.

Technical Specifications

- Height 156mm
- Width 651mm
- Depth 180mm

Product Code:

Dispenser: ACTDR01
Cover Roll: ACTR01

Features

- Suitable for Changing Places applications
- Easy 6 screw fixing
- Sleek stainless steel appearance

Need more information?

 web email call

CHS Wipe Dispenser



Accessories

The CHS wipe dispenser is a wall mounted, super sleek, streamlined and fully functional in various applications. The stainless steel wall mounted stand works with many equipment wipes and is easy to refill. The unique sturdy Silicon valve allows a single wipe to be dispensed minimising wastage as wipes are dispensed.

Technical Specifications

- Overall Height - 300mm
- Overall Width - 240mm
- Overall Depth - 230mm

Product Code: CHSWWD01

Features

- Easy fastening
- Modern and sanitary stainless steel
- Wipes not included

Need more information?

web

email

call

CHS Shower Seat

Accessories



The CHS shower seat is a stainless steel construction and a plastic top. This shower seat has been specifically designed to meet all Changing Places seat requirements.

Technical Specifications

- Overall Height - 480mm
- Overall Width - 960mm
- Overall Depth - 450mm
- 960mm long seat
- 450mm depth to the wall (40mm max gap permitted)
- Slip resistant
- Front corners that are rounded to a radius of 10-15m
- Top edges that are rounded with a minimum radius of 2-3mm
- Drainage holes where provided between 4-6mm
- Fastenings, materials and construction of the seat shall be able to withstand a force of 1100N (110Kg) in any direction at any position and in any direction without failing or loosening of fastenings.





Product Code: CHSSS450D

Need more information?

web

email

call

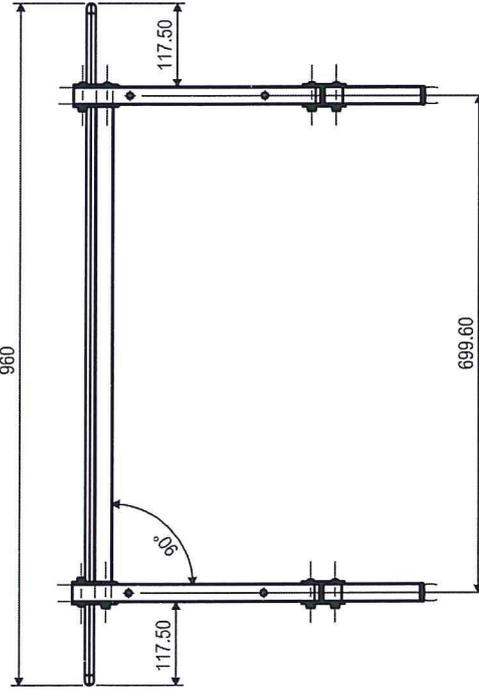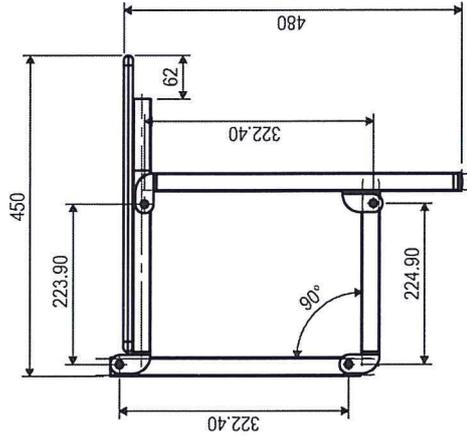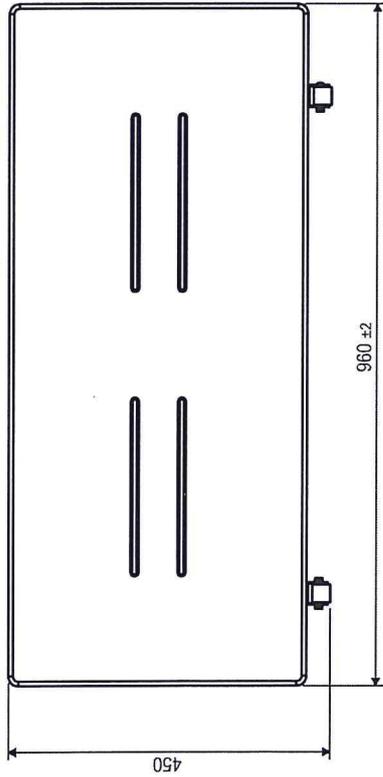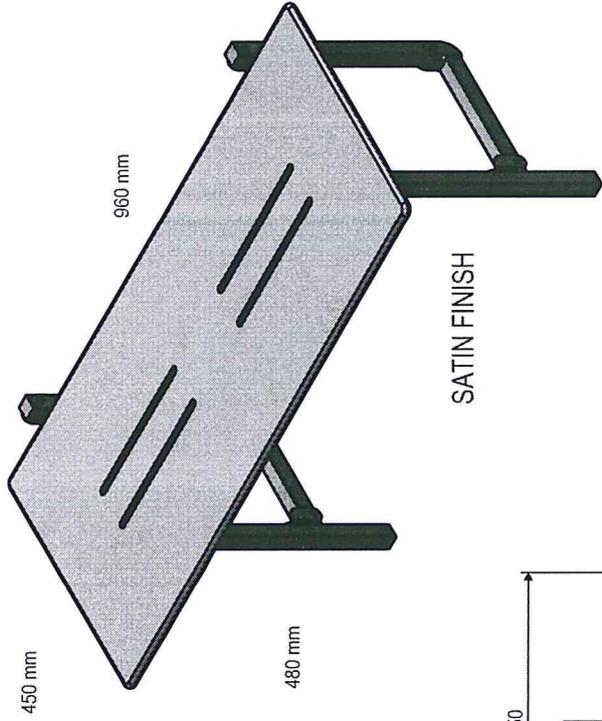

		REV	0	DATE	APR032017	Changes	ORIGINAL DRAWING	ECN #		Dwn	NILO	Chk	DON	SHOWER SEAT 960 x 450 mm <small>UNLESS OTHERWISE SPECIFIED: DIMENSIONS ARE IN MILLIMETERS TOLERANCES: LINEAR: ± 1.0mm; ANGULAR: ± 0.5°</small>		REF	
		<small>MATERIAL ALTERNATIVES: FINISH: WEIGHT:</small>															
<small>A4 SCALE 1:10 SHEET 1 OF 2</small>																	

Ropimex RKL

Curtain Rails



With this space-saving model, you can create a screen system in seconds that can be used for washing, changing clothes, as well as for radiation or inhalation treatment, as a corner or U-shape solution.

Technical Specifications



Features

- Can be assembled in seconds
- Compatible with Ceiling Hoist System unlike many other rail and curtain systems
- Simple and easy to use with many varying installation features to suit every situation
- Optional curtain styles and colours to match or brighten any room

Corner solution

Extension

Telescopic screen cubicle

max. 120/115 cm
min. 75/ 70 cm

Telescopic screen cubicle

max. 150/145 cm
min. 90/ 85 cm

Need more information?

web

email

call

Ropimex Curtains – TCS Trevira Curtain Rails



Ropimex® curtains made from Trevira®CS / Trevira®CS BIOACTIVE
Our flame-retardant curtains with a feel-good factor. For that extra safety.

Applications:
As curtains to protect privacy in healthcare and care facilities, and as attractive room partitions
As a room divider covering

TCS ropimex® curtains made of Trevira®CS have been tried and tested for 20 years



Trevira®CS has passed every fire test worldwide and fulfils the fire test standards of all the important industrial nations.

The woven-in company logo ropimex® with the Trevira®CS label guarantee you this safety for ropimex® curtains.

Trevira®CS has an antibacterial effect. The physiological and dermatological properties are excellent.

Trevira®CS Bioactiv is also **permanently** antimicrobial.

Curtains made of Trevira®CS are washable up to 60° C with disinfectant.

Size table

	Width	x	Height in cm
TCS 110	110	x	175
TCS 121	140	x	120
TCS 122	140	x	140
TCS 123	140	x	175
TCS 151	170	x	120
TCS 152	170	x	140
TCS 153	170	x	175
TCS 200	200	x	175
TCS 211	230	x	120
TCS 212	230	x	140
TCS 213	230	x	175
TCS 215	250	x	175
TCS 217	270	x	175
TCS 300	300	x	175
TCS 330	330	x	175
TCS 360	360	x	175
TCS 390	390	x	175
TCS 420	420	x	175
TCS 460	460	x	175
TCS 500	500	x	175

Available in 11 colours



lilac



heather



sun yellow



melon



peach



azure



light blue



white



light green



rose



champagne

Trevira®CS Bioactive available in 4 colours: white, light green, light blue and peach.

Need more information?

web

email

call

Ropimex RFW



Curtain Rails

Can be used flexibly as a screen or partition in a hospital or care facility, this robust and impact-resistant folding wall can be brought into position swiftly and just as quickly returned to its storage position on the wall. The number of leaf sections can be extended later as required.

Technical Specifications

- 3 Heights from 145 – 185cm
- 3 to 10 leafs possible
- Up to 360° rotation

Features

- Space-saving, fold-away and ready to use again in seconds
- Non-see through but translucent material create a pleasant and light atmosphere
- Easy clean and scratch resistant surfaces
- Modern designs and colours to add an individual touch

Need more information?

 web email call

Ceiling Hoist Help Guide



**Changing Places
Transforming Lives**

1 Product description:

This ceiling hoist is intended for lifting children and adults up to 200kg. This equipment must not be used for any other purpose.

2 Pre use checks:

Hold and pull the carry bar, moving the hoist to an area free of obstructions. Take care to familiarise yourself with the functions and controls (see fig 1).

3 Instructions for use

The hoist is controlled via the hand control (see fig 2). Do not pull the hand control with force. When the hand control is not in use, clip over the carry bar. Keep all obstructions clear during lifting.



4 Slings

This hoist is designed for use with loop style slings only. If unsure do not use your sling.

5 Trouble shoot

If the hoist is not working check to see if the emergency stop cord has been pulled into the down position, disabling the hoist (see fig 3).

To re-enable push white clip, attached to the top of the red cord, up toward the hoist. If unable to reach use the carry bar.

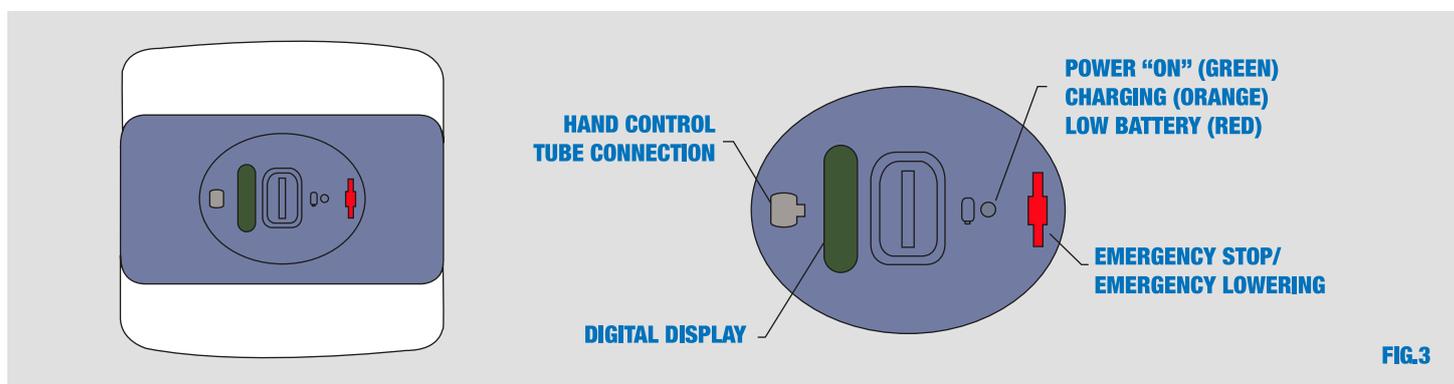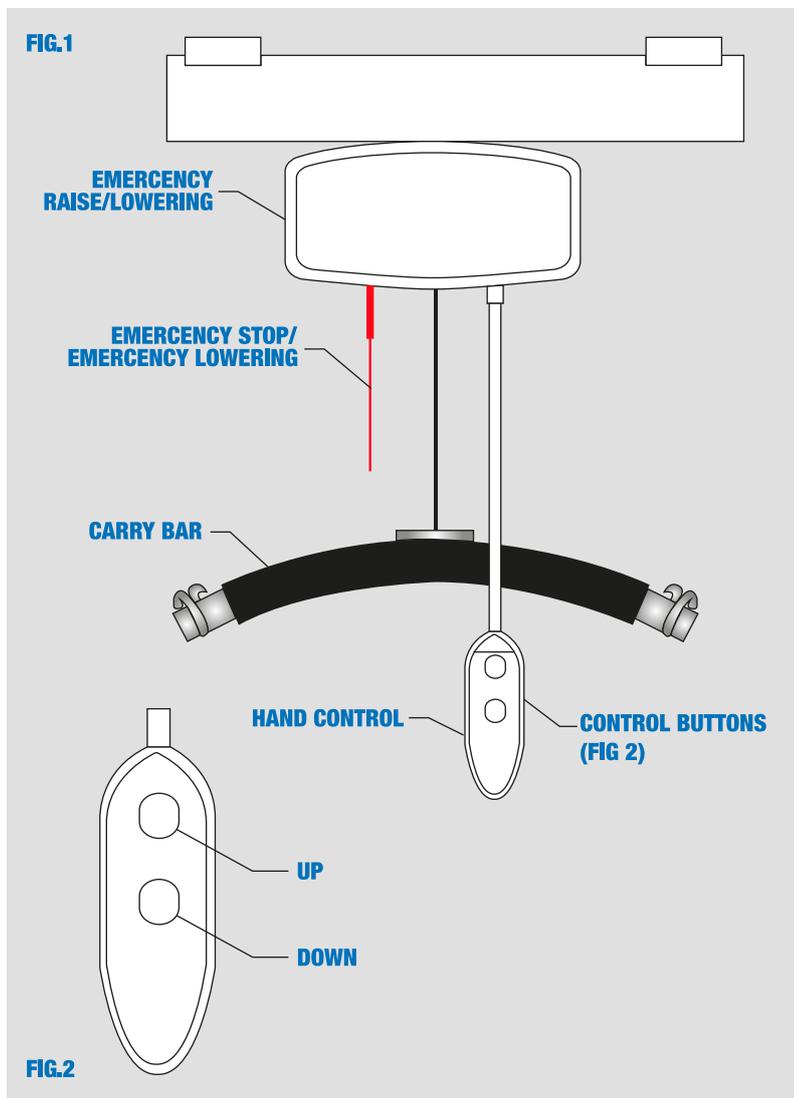

6 Emergency lowering

If the hoist becomes stuck or won't lower, pull the red cord and press the down button on the hand control simultaneously to lower the carry bar to the desired height. Use this function only in an emergency.

Ceiling Hoist Help Guide



Changing Places
Transforming Lives

1 The STELLAR Ceiling Hoist. What to do & how it works?

- Hold and use the spreader bar to pull the hoist into a clear area for safe use. Make sure there are NO obstructions around you.
- Take care to familiarise yourself with the features of this hoist.
- Max user weight 200kg.

2 The hand control, spreader bar and your sling

- The hoist has 2 buttons - up and down
- Do not pull hard on the hard control cord
- When the hand control is not in use clip the hand control over the spreader bar (left) to keep it safely out of the way
- Keep all obstructions clear when lifting

3 Emergency STOP & Lower

- On the hoist is **red emergency** stop cord. This cord is only for emergencies. Pulling this cord will stop power to the hoist and the hoist motor will stop sharply.
- To **emergency lower** pull and keep hold of the chord at the lowered position.
- To **reset** the switch press the metal button back up.

4 Slings

- Unfasten slings to disengage after emergency stop & lowering.



Need assistance? Call CHS Healthcare

1300 789 420

www.chshealthcare.com.au

Proudly supported by



CHS Healthcare
We're Here For YOU

Change Table Help Guide



**Changing Places
Transforming Lives**

1 Product description:

This changing table is intended for showering and changing children and adults weighing up to 200kg, the table must not be used for any other purpose.

2 Pre use checks:

Check for obstructions or items that may become trapped when the changing table is operated. Check the changing table will raise and lower properly using the handset.

3 Lowering the bed:

To fold down the bed to the usable position, simply hold the bed and lower as indicated by the arrow in Fig 2. To stow bed, lift and fold bed upwards until parallel with cover.

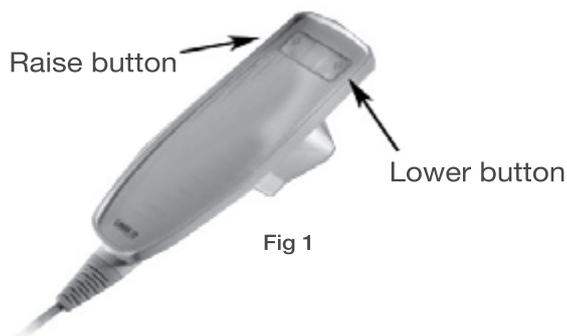



Fig 1




Fig 2

4 Instruction for use:


Operation of the change bed is via the two button hand control. Position the bed at a comfortable height. See Fig 1

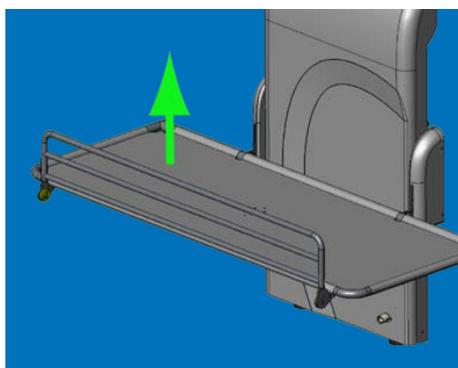

Operating the guard (where fitted). Folding the guard. Lift the guard upwards.

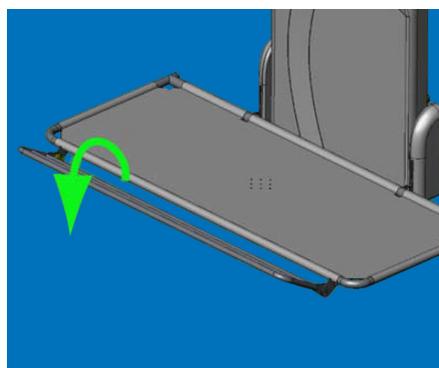

Then rotate it towards you and then fold under the bed.

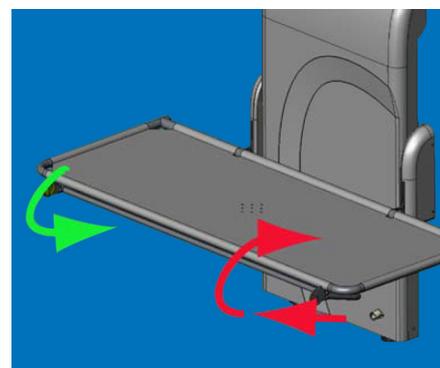

It will then lock under the bed. To replace the guard lift it slightly towards you, it will then release, rotate upright and drop into position.

Need assistance? Call CHS Healthcare
1300 789 420
www.chshealthcare.com.au

Proudly supported by  **CHS Healthcare**
We're Here For You

Change Table Help Guide



**Changing Places
Transforming Lives**

1 Product description:

This changing table is intended for showering and changing children and adults weighing up to 200kg, the table must not be used for any other purpose.

2 Pre use checks:

Check for obstructions or items that may become trapped when the changing table is operated. Check the changing table will raise and lower properly using the handset

3 Lowering the bed:

To fold the bed to the useable position, simply hold the bed and lower as indicated in figure 1. To stow the bed lift and fold bed until parallel with the cover.

4 Instruction for use:

Operation of the change bed is via the hand control. See figure 2.
To use the guard, pull on it from the stowed position and rotate it around to the upright position. Once upright it will drop and lock in place (see figure 3). To stow, lit up until it can be rotated backward. Rotate under the bed until the magnets hold it.

